




HALL EFFECT APPARATUS



SKU: frm-hall-ge

Category: [Scientific Equipment](#)

GALLERY IMAGES

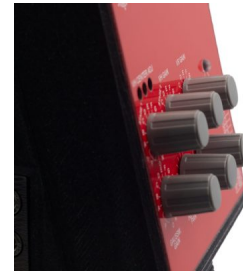

PRODUCT DESCRIPTION

Measurements:

- Resistance vs. T
- Hall voltage V_H vs. T and B
- Hall Mobility μ_H vs. T
- Evaluation of dopant density N_a
- Energy gap E_g^0 linearly extrapolated to T=0K
- Mobility ratio $b=\mu_e/\mu_h$
- Earth magnetic north pole

Introduction:

The Hall effect is important in the investigation of the properties of the materials because it allows direct measurements of the free carriers concentration (with their sign) both in metals and in semiconductors.

It is also technologically important because it is the working principle of modern magnetic field sensors.

The Hall effect essentially is due to the Lorentz force F acting on each electric charge q moving with velocity v in a magnetic field B .

While in principle the experiment is not overly complex, in practice it often is. This apparatus solves in advance a great number of challenges that would impair the correct experiment and take weeks if not months to sort out:

- Sourcing, cutting and polishing a sample with the appropriate parameters to conduct the experiment successfully
- Stable ohmic electrical contacts to the sample
- Isothermal environment for the sample avoiding thermal gradients
- Controlling the temperature in a wide range -150°C to +150°C
- Providing a stable and precise constant-current to the sample
- Providing a variable, stable, and uniform magnetic field to the sample
- Keeping self-heating due to the bias current to a minimum
- Providing adequate amplification to all the signals with a variable/programmable gain
- Measuring the magnetic field
- Keeping in check errors:

- Microphonic pickup of the sample wiring inside the magnetic field
- Uneven magnetic field
- Uneven thermal gradients
- Seebeck effect
- Parasitic Schottky diodes on the samples and the likes
- Amplifier drift, offset
- Unbalanced bias currents of the differential amplifiers
- Reference voltage stability and precision
- Constant current generator drift and precision
- ...
- Providing a mean to compensate for imperfections in the sample contact positioning that would overwhelm the delicate V_h signal
- Measure the resistance of the sample through a 4-wire measurement system
- Measuring the sample temperature precisely (not easy at cryogenic temperatures)
- Shutting down the heater when required to avoid damage
- Safe power supply
- Providing adequate shielding to the measurement to avoid noise pickup

This apparatus is a 5th generation model:

- 1st and 2nd generation from 1990 to 2015
- 3rd from 2015 to 2017
- 4th from 2017 to 2019
- 5th 2020+

Specifications:

Electronics:

- Compact and all-in-one
- 5-point hall voltage measurement with virtual centre point balance (to compensate for the sample contact positioning)
 - balanced Hall voltage output, unbalanced Hall voltage output 1, unbalanced Hall voltage output 2
- High precision amplifier stages for V_h , V_r :
 - 10uV offset voltage

- 1 to 4096 programmable gain with 0.1% precision
- 1.5nA input bias current
- Gaussmeter probe with 5% precision, variable gain from 1 to 4096, stable 1mA bias current supply
- Two power-level heating system
- Stable 0-20mA constant bias current generator with 0.1% readout on a 0-5V output
- Lock-in modulation input
- 12V 45W power supply
- Grounding plug/crocodile clip
- K-type thermocouple with cold-point compensation circuitry
- Safety 150°C heater cut-out with 100°C alarm reset
- Safety alarm at 60°C "HOT - do not touch!" with blinking LEDs
- Retention of the configured parameters

Mechanics:

- Variable magnetic field
 - Permanent rare-earth magnets
 - Up to 0.5T field
 - soft-iron toroid
- -150°C to 150°C with liquid nitrogen and internal heater
- Anodized aluminium base and supporting rods

Accessories:

- Gaussmeter probe
 - 1mA bias current
 - up to 5T measurement
- Gaussmeter positioning jig: allows the positioning of the gaussmeter probe at the center of the sample

Included in the kit:

- Aluminium base and holder
- Permanent magnets and soft iron toroid with a variable magnetic gap
- Isothermal sample holder with heater, thermocouple, cold finger, Ge-P sample
- Stainless steel dewar (for liquid nitrogen or acetone+cold ice)

- Advanced electronics with variable gain, high-impedance, low noise, low drift, low bias current differential amplifiers; constant current generator, gaussmeter control circuitry, heater control circuitry, thermocouple circuitry
- Gaussmeter probe and connection wire
- Gaussmeter positioning jig
- Two DE-15 to 5-BNC cables
- 12V 45W power supply
- Digital Handbook

Needed (not included in the kit) :

- Datalogger (recommended Interface LabPro or equivalent)

Innovative Handbook:

- [PDF](#)
- [HTML](#) for projection in the classroom with variable font size, background colour.
- [Open source](#): The professor/student can contribute feedback and modifications to the handbook
- Included interactive data analysis with open-source R language
- Constantly updated